

**DST-6**

**Desktop Electronic Counter**

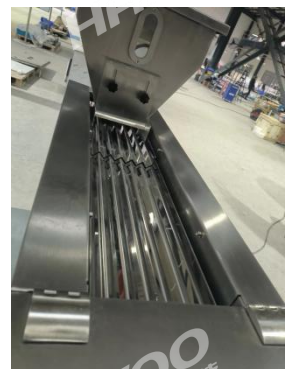
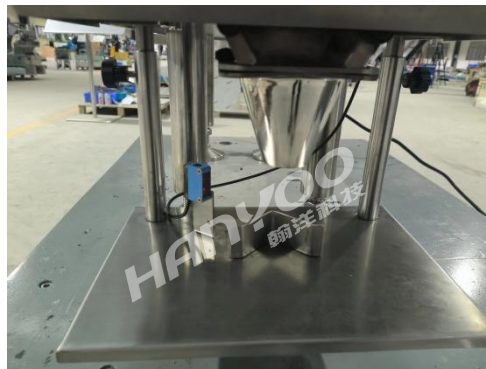


**TECHNICAL  
PARAMETER**

**南京翰洋智能设备有限公司**  
NANJING HANYOO PHARMATECH CO., LTD.



# DST-6 Desktop Electronic Counter



---

## INTRODUCTION

The empty bottle is placed in the bottling position manually, and the drug enters the drug bin in an orderly manner through the vibration of the feeding corrugated plate. The drug bin is equipped with a counting photoelectric sensor, and the drug falling into the drug bin is quantitatively counted by the photoelectric counting sensor and then loaded into the bottle. Manually filled bottles are removed and new empty bottles are placed to be filled.

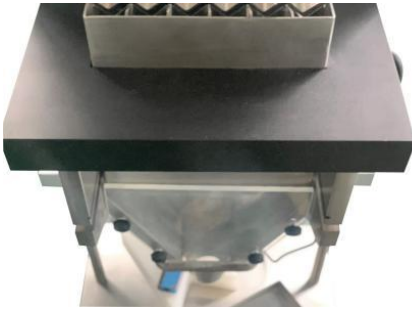
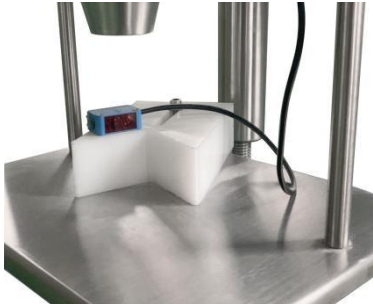

## PRODUCT FEATURE

1. Strong compatibility, solid preparations or solid particles such as tablets (including shaped tablets), capsules, soft capsules (transparent, opaque), pills can be counted and bottled;
2. Strong applicability, the use of the company's anti-high dust photoelectric sensing technology, in the case of high dust can work stably;
3. Vibrating type feeding, unique flip plate subassembly mechanism feeding, discharging uniform and stable, no damage to drugs;
4. The leak nozzle of special structure can prevent blocking of medicine, and the medicine bottle with small mouth can also be quickly filled (optional);
5. High degree of intelligence, with a number of detection alarm control functions such as no bottles, abnormal medicine bottles or automatic elimination of certain size of broken particles, fault self-test;
6. Easy to replace varieties, disassembly and clearance without any tools, convenient and fast;
7. Touch screen control panel, easy to adjust and use;

## TECHNICAL PARAMETER

MODEL	DST-6
Scope Of Application	$\Phi$ Pill $\leq 20\text{mm}$ ; $5.5 \leq \Phi$ Tablet $\leq 22\text{mm}$ ; Capsules; Softgel
Suitable Bottle Type	Bottle diameter 20-100mm, height 40-180mm
Hopper Capacity Type	About 11L
Production Rate	10-30 bottles/min
Power	220V/50HZ 0.6KW
Size	1090*280*840
Weight	100KG

## DETAIL INTRODUCTION

	The materials pass through the vibrating plate and fall into the medicine bin by the counting photoelectric system
	The single discharge port bottle inspection photoelectric eye can determine whether there are bottles at the discharge port (it won't work without bottles).
	The materials move orderly into the medicine warehouse through the vibration of the corrugated plates

## MAIN CONFIGURATION LIST

NO.	NAME	MANUFACTURER
1	Touch screen	Delta-China
2	Photoelectric sensor	Panasonic
3	Count the Photoelectric Head	Adopt Everlight Electronics., Ltd
4	Electric leakage protection	Schneider (France)
5	main control board	Imported high-speed chips
6	Switching power supply	Schneider (France)



[Http://www.hanyoopm.com](http://www.hanyoopm.com)